

Product: Chicken Anti-GFP protein, affinity purified.

Catalog Number: GFP

Lot Number: 70102

Size: 0.1ml, 0.5mg/ml in PBS with 0.02% Na Azide

Immunogen: Purified GFP, fusion protein.


Specification: The antibodies bind specifically to the GFP protein in cell lysate containing GFP fusion protein in a western blot analysis.

Suggested Usage: The product can be used at 1:2,000 to 1:5,000 for both immunoblot and ELISA testing. User needs to try out the titrations for optimization in immunocytochemistry (ICC) and immuno-histochemistry (IHC) assay. The product can be used in combination with GeneTel's anti-IgY conjugates such as RCYHRP (anti-IgY HRP) and RCYAP (anti-IgY AP) for ELISA and Western Blot analysis, and RCYRP or RCYFITC for ICC/IHC assay.

Short Term Storage: 2-8°C for one week.

Long Term Storage: -20°C as an undiluted liquid for up to two year.

Please feel free to call us if you have any question regarding to this product or any other GeneTel's product.

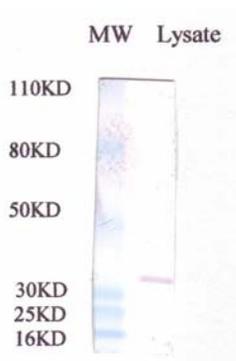

Western Blot analysis of cell lysate containing GFP protein with chicken anti-GFP at a 1:4,000 dilution.